NOTE: The following information is in addition to the owner's manual supplied with your scooter. It explains how to read your battery voltage and the approximate remaining charge of the batteries using the voltmeter and how to perform basic troubleshooting using the diagnostic codes. You must read the scooter owner's manual and all supplemental information before use. Retain this information for your records, and contact your authorized Pride Provider with any questions.

Reading Your Battery Voltage

The battery condition meter on the tiller console indicates the approximate strength of your batteries using a color code. From right to left on the meter, green indicates fully charged batteries, yellow indicates a draining charge, and red indicates that an immediate recharge is necessary. See figure 1. To ensure the highest accuracy, the battery condition meter should be checked while operating your Scooter at full speed on a dry, level surface.

HOW TO READ LITHIUM BATTERY USAGE ON YOUR VOLTMETER:

- 3 Green LEDs = 90% - 100% Charge
- 2 Green LEDs = 80% - 90% Charge
- 1 Green LED = 60% - 80% Charge
- 4 Yellow LEDs = 50% - 60% Charge
- 3 Yellow LEDs = 40% - 50% Charge
- 2 Yellow LEDs < 30% Charge
- 1 Yellow LED < 20% Charge
- All Red LEDs = Dead Battery

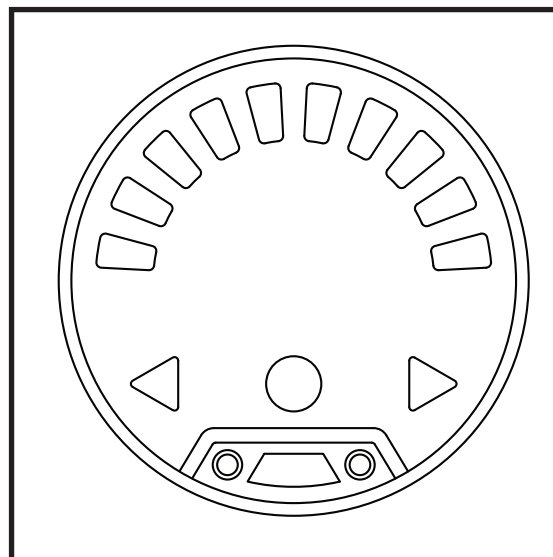
Figure 1. Battery Condition Meter

DIAGNOSTIC CODES

The diagnostic codes for your Scooter are designed to help you perform basic troubleshooting quickly and easily. A diagnostic code will sound in the event one of the conditions listed below develops. A diagnostic code will identify the condition will repeat four (4) times and then stop beeping. To repeat the diagnostic code, remove the key, then reinsert it and turn it to the “on” position.

NOTE: *Your Scooter will not run unless the diagnostic code condition is resolved and the Scooter has been turned off, then turned back on.*

DIAGNOSTIC CODE	CONDITION	SOLUTION
■ (1) (steady tone)	Battery voltage is too low to operate the Scooter.	Charge batteries as soon as possible.
■■ (2)	Battery charge is too low.	Charge batteries.
■■■ (3)	The Scooter’s battery voltage is too high to operate the scooter, or the charger is still connected to the off-board charger port.	Unplug charger and/or turn key off, then back on again. If traveling down a slope, reduce Scooter speed to the minimum setting.
■■■■ (4)	Hall effect sensor/Motor Error	Turn the Scooter off for a few minutes, then turn your Scooter back on.
■■■■■ (5)	The manual freewheel lever is in the (forward) freewheel position.	Remove the key from the key switch, then push the manual freewheel lever to the drive (rearward) position and restart your Scooter.
■■■■■■ (6)	Throttle control lever not at center position at start up.	Verify that the batteries are fully charged. If the batteries are fully charged, then return the throttle control lever to center position, turn Scooter off then back on.
■■■■■■■ (8)	Motor volts error.	Contact your authorized Provider for assistance.
■■■■■■■■ (9)	Other internal errors.	Contact your authorized Provider for assistance.

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