

www.coretechortho.com

coretech[®]
by **vive health**


OWNER'S MANUAL
906 MOTION ANKLE BRACE
SKU: SUP2044



To see all FAQs in one place, visit, [coretechortho.link/17h](https://www.coretechortho.com/link/17h)

906 MOTION ANKLE BRACE

1. Loosen and open all straps.

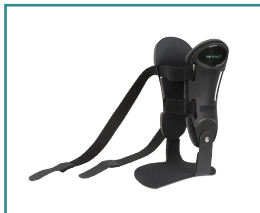

2. Insert the plastic base into your shoe.



3. Insert your foot through the motion ankle and into your shoe.



4. Secure the bottom strap around your calf first, then the top strap.



Test for comfort and fit, adjust as necessary.



CARE INSTRUCTIONS

The side pads are removable for cleaning. Remove them and hand-wash in cold water with mild soap, rinse, and air dry. Wipe the plastic with a damp cloth and dry immediately. Do NOT machine wash or dry. Do NOT bleach.





Distributed by

vive
health

8955 Fontana Del Sol Way
Naples, FL 34109
1-800-487-3808
www.vivehealth.com