

www.coretechortho.com

coretech[®]
by vive health

OWNER'S MANUAL
639 LSO
SKU: SUP3016BLK



To see all FAQs in one place, visit, [coretechortho.link/cmw](https://www.coretechortho.com/cmw)

SUP3016BLK

639 LSO

WHAT'S INCLUDED

- 1x 639 LSO
- 1x Top Panel
- 1x Bottom Panel

CARE INSTRUCTIONS

Spot clean with a damp cloth, air dry. Do not machine wash or dry.



WARNINGS

The product is intended for application by a licensed health care practitioner as directed by a physician or other qualified medical authority. Permanent injury may result from use other than as directed by a physician or other medical authority.

ADJUSTING THE BRACE

Before putting the brace on, it's important to adjust the brace to a suitable size. The upper left and upper right belts are labeled so that the user can identify the correct orientation when adjusting. There is a size range representing waist size printed on the inside of the belt.



BELT SIZING:

1. Disconnect the fabric fasteners on each side of the belt from its tab on the back panel.
2. Fold each side of the belt to the desired size. The fold in the strap will determine the functional length of the straps.



****You may cut the belt to the desired length as well.****

3. Reattach the straps to the fastener tabs on the back panel.




ANTERIOR SUPPORT SIZING:

1. Open the fabric cover of the bottom anterior panel and loosen both screws.



2. Place the anterior support on the patient with the bottom panel around the abdomen and the top panel in between the middle to bottom of the sternum.

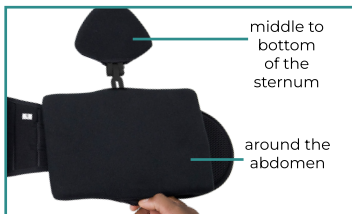

3. Note the length and tighten the screws to secure in place.

SUP3016BLK

4. You can also adjust the angle of the top anterior panel for further support. Loosen the two screws located on the side of the hinge.
5. Place the support on the patient and adjust to the desired angle. Tighten the two screws back in place to secure.



APPLICATION

1. Position the brace so that the back panel is right side up and centered on your spine, just above the waist and the size labels on the inside are facing right side up.

****Assistance May Be Required****



2. While keeping the back panel centered and both belt straps open, position the anterior support to the patients front and wrap the belt to secure in place.



3. Grab the compression cords using the thumb loops and slowly pull the cords in unison away from your body until the desired amount of compression is achieved. Secure to the front of the brace.



4. Adjust the support as needed for the most comfortable fit.



Distributed by

vive
health

8955 Fontana Del Sol Way
Naples, FL 34109
1-800-487-3808
www.vivehealth.com