

www.coretechortho.com

coretech[®]
by vivehealth

OWNER'S MANUAL
ROM WALKER BOOT - TALL
SKU: SUP3037BLK



To see all FAQs in one place, visit, coretechortho.link/aib

ROM WALKER BOOT - TALL

PUTTING THE WALKER ON



For a video demonstration,
check out coretechortho.link/gbt

1. Remove the velcro covers on each side of the boot before use to secure the liner in place.



2. Open the inner liner and insert your foot. Secure the bottom then the top to a comfortable fit.



3. Secure the outer straps from bottom to top and tighten to a comfortable fit. Do not overtighten.



Optional: If the boot is uncomfortable, loosen the straps and use the extra pads included to pad the inner liner or straps.




4. Set the air valve to "INFLATE" by turning it clockwise.



SUP3037BLK

5. Press the pump by hand to inflate to your desired compression.
6. Push the tabs on the ROM dial to set your flexion and extension angle.
7. When removing the boot, set the air valve to "DEFLATE" by turning it counter-clockwise. Use hand compression to speed up the deflating process if needed.

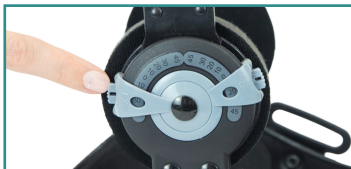

CARE INSTRUCTIONS

- Remove the padded liner. Wipe the booth with a damp cloth and dry immediately.
- Do NOT machine wash or dry

coretech[®]
by **vive health**

Distributed by

vive
health

8955 Fontana Del Sol Way
Naples, FL 34109
1-800-487-3808
www.vivehealth.com