



CRUTCH PADS WITH ERGONOMIC POUCH

CSH1104



HAVE MORE QUESTIONS?

Check out our list of Frequently Asked Questions at
vivehealth.link/klm for helpful answers.



And if that doesn't answer your question, our
customer service team would love to help!
Feel free to connect with them by phone,
e-mail, or chat on our website

service@vivehealth.com
1-800-487-3808
vivehealth.com



Distributed by



110425



To see all the FAQ's
in one place,
visit vivehealth.link/klm

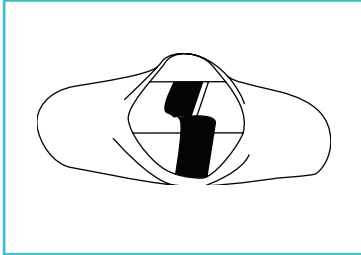
OWNER'S MANUAL

CRUTCH PADS WITH ERGONOMIC POUCH

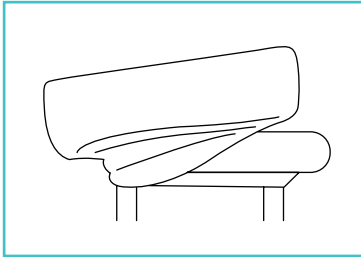
QUICK START GUIDE

Step 1:

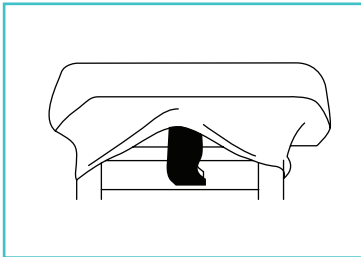
Open up the fastening straps.

**Step 2:**

Stretch cover over the crutches.

**Step 3:**

Attach the fastening straps for a snug fit.



POUCH INSTALLATION

Follow the same steps to install the extended pouch. Use the zipper and the Snaps to secure as needed.

MAINTENANCE & CARE

- Remove cloth material from foam padding. Machine wash in cold water on a gentle cycle.
- Use mild detergent.
- Do not use bleach.
- Tumble dry on low heat or hang to air dry.
- Do not iron.