

vive
vivehealth.com

Distributed by

vive
health

8955 Fontana Del Sol Way
Naples, FL 34109
1-800-487-3808
www.vivehealth.com

FOREARM CRUTCHES

MOB1014




Check out this how-to video
at vhealth.link/3cd

vive
vivehealth.com

QUICK START GUIDE

FIND THE CORRECT HEIGHT

1. Let your arm hang naturally by your waist. Adjust the forearm crutch height so the handle touches your wrist.

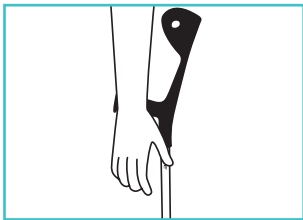

CANE HEIGHT ADJUSTMENT

1. Press in the button snap.
2. Slide extension tube up or down to desired height. Button snap should secure cane into place.



WARNINGS ⚠

- Make sure the brass pin is locked into place before use.
- Weight Capacity: 300 LBS