

LED FOLDING CANE

MOB1017



To see all the FAQs,
visit [vhealth.link/64dd6](https://vivehealth.link/64dd6)

vive
vivehealth.com

QUICK START GUIDE


QUICKSTART

1. Unbutton the wrist strap that is holding the cane together.
2. Hold the handle and release the folded sections. The cane will automatically snap into place.
3. Untwist the LED light cap and remove the battery protection tab to use the LED light (see changing the batteries section for help).



CANE HEIGHT ADJUSTMENT

1. Loosen the locking latch.
2. Adjust brass pin to desired position.
3. Tighten the locking latch.

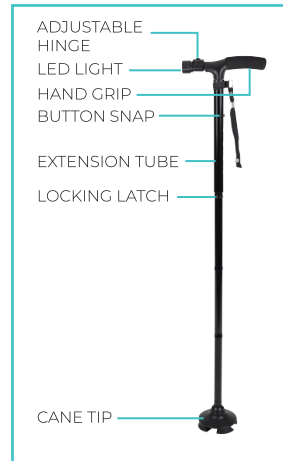



Distributed by



8955 Fontana Del Sol Way
Naples, FL 34109
1-800-487-3808
www.vivehealth.com

FIND YOUR CORRECT CANE HEIGHT

Let your arm hang naturally by your waist. Adjust the cane height so the handle touches your wrist.

CHANGING BATTERIES

1. Twist cap counterclockwise and remove the old batteries.
2. Insert 4 new LR44 Batteries negative side first into the black casing. Place the black casing back into the battery holding slot, negative side first.
3. Twist cap clockwise to lock into place.



Note: Twist the cap all the way on to ensure the batteries are engaged.

WARNINGS ⚠

- Button batteries can be dangerous and cause serious injuries, including death, if swallowed or inserted into a child's nose. Button Batteries can cause potential choking hazards. Please keep away from small children.
- Make sure the brass pin is locked into place before use.
- **Weight Capacity:** 250 LBS