

DUAL TRIGGER FINGER SPLINT

SUP2091BLK

vive
vivehealth.com



To see all FAQs in
one place,
visit vhealth.link/hks

OWNER'S MANUAL

DUAL TRIGGER FINGER SPLINT


WHAT'S INCLUDED

- Dual Trigger Finger Splint



APPLICATION

1. Open up the splint and insert your thumb into the opening.

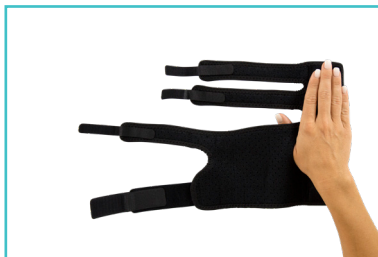

2. Secure the four straps starting at the wrist.

NOTE:

There are different ways to secure fingers depending on your needs.



3. To use the splint on the left hand, open the flap covering and remove the splint, then insert it into the opposite side.



4. Once the metal splint has been moved to the other side, adjust the hook and loop Y-Straps to use on the left hand.



5. Once all 4 straps and the splint have been moved to the opposite side apply steps 1 and 2 on the left hand.

NOTE:

The Bottom Strap that secures around the wrist is trimmable if too long. Just remove the Y-Strap and cut to desired length.

CARE

- Remove splint, then hand wash with cold water and a mild detergent.
- Hang dry.

WARNINGS

- Consult your doctor before use.
- Stop if pain worsens.
- If you experience any pain, swelling, sensation changes, or any unusual reaction while using this product please consult your doctor immediately.
- Assistance may be required during the application of this support.
- Do not over tighten, doing so may lead restricted blood circulation.

NOTES



HAVE MORE QUESTIONS?

Check out our list of Frequently Asked Questions at
[vivehealth.link/hks](https://vivehealth.com/healthlink/hks) for helpful answers.



And if that doesn't answer your question, our
customer service team would love to help!
Feel free to connect with them by phone,
e-mail, or chat on our website

service@vivehealth.com
1-800-487-3808
vivehealth.com

Distributed by



8955 Fontana Del Sol Way
Naples, FL 34109
Made in China